




# MCS Micro-cleaning System

The MCS is a micro-valve plug-in system for the MM3A-LMP micromanipulator. It enhances your system by allowing non-destructive removal of unwanted particles from sensitive surfaces.

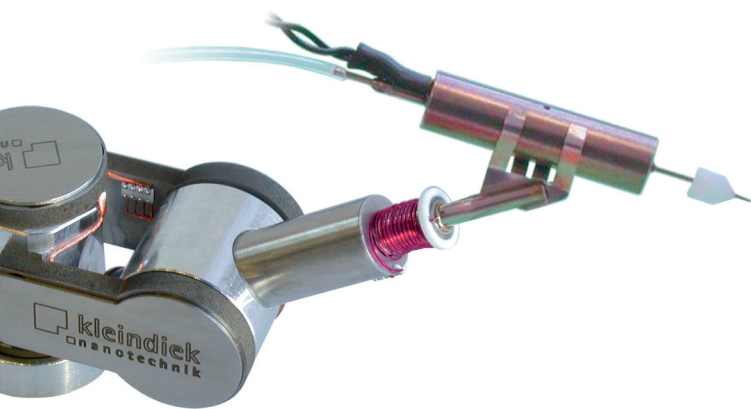

The unique combination of precision manipulation and controlled application of compressed air provides you with a reliable, fast and easy-to-use cleaning system. Control over the strength and length of each pressure impulse is possible.

The plug-in tool can be quickly and easily retro-fitted to your existing MM3A-LMP micromanipulator.

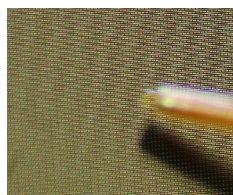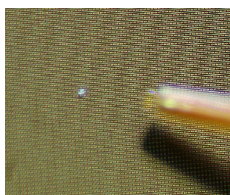

## Applications

- Micro-cleaning

## Technical specifications

- Microvalve length 31.4 mm
- Microvalve diameter 5.6 mm
- Inlet diameter 0.6 mm
- Outlet outer diameter 0.5 mm
- Outlet inner diameter 0.25 mm
- Weight 3 g

## Further information

- Contact us at [info@kleindiek.com](mailto:info@kleindiek.com)
- Find your local agent at [www.kleindiek.com](http://www.kleindiek.com)

