



# MM3A-EM Micromanipulator

Over the last few years, the MM3A-EM micromanipulator has created an unparalleled new dimension of quality in the field of Electron Microscopy. It is employed in a wide spectrum of SEM, FIB and other microscopes for an even wider range of applications and it has become the industry standard for OEM and retrofit solutions.

**Give your microscope a hand:** use the MM3A-EM to add new capabilities and functionality to your instrument.

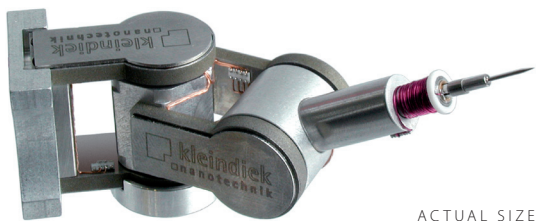

ACTUAL SIZE

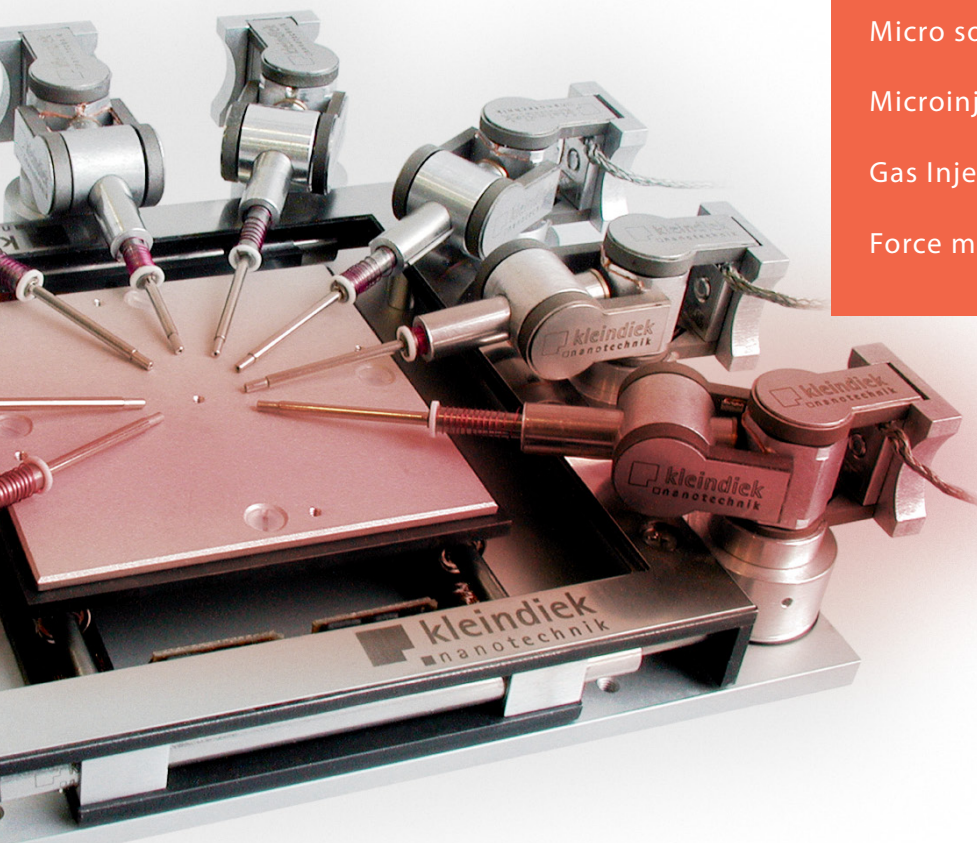

## APPLICATIONS

- Electrical probing (FA)
- Nanomanipulation
- In-situ lift-out
- Materials science

## PLUG-IN TOOLS

- Micro four-point probe
- Microgripper
- Rotational tip
- Low current measurement kit
- Micro soldering unit
- Microinjector
- Gas Injection System
- Force measurement system

# MM3A-EM Micromanipulator

## More compact and more flexible

- Small and practical
- Plug-and-play system with modular components
- Interfacing solutions for most SEM/FIB instruments (including load lock)
- Fast setup and removal
- Effortless work with multiple manipulators
- Useful plug-in tools

## Clearer and simpler

- Result-oriented operation which leads to increased throughput
- Intuitive control interfaces and software
- User-friendly and easy to learn
- Quick and easy probe tip exchange
- Compact, stand-alone electronics
- Pioneering cabling technology with compact vacuum feedthrough

## More robust and more stable

- Compact construction delivers higher resonance frequencies
- Excellent stability
- Low drift (1 nm/min)
- Reliable operation (one year endurance test)
- Virtually unsusceptible to vibrations
- Fast pre-positioning by hand

## Faster and more precise

- No backlash or reversal play
- Sub-nanometer resolution (0.25 nm)
- Extensive working range (100 cm<sup>3</sup>)
- No "blind axis" like with cartesian systems
- Coarse and fine displacement in one drive
- High operating velocity (up to 10 mm/sec)

## Technical specifications

- Length 62.1 mm
- Width 20.4 mm
- Height 25.4 mm
- Weight 45 g
- Operating range AB 240°  
Operating range C 12 mm
- Speed AB up to 10 mm/s  
Speed C up to 2 mm/s
- Resolution A 10<sup>-7</sup> rad (5 nm)  
Resolution B 10<sup>-7</sup> rad (3.5 nm)  
Resolution C < 0.5 nm
- Holding force 1 N
- Holding torque 3 to 4 Nmm
- Lift Y 5 g
- Probing current range 10 nA to 100 mA
- Maximum probing voltage 100 V
- Probing signal resistance 7.0 Ω
- Temperature range 273 K to 353 K  
UHV version 273 K to 393 K
- Lowest pressure 10<sup>-7</sup> mbar  
UHV version 2 × 10<sup>-10</sup> mbar
- Mounting M4 tapped hole
- Material Stainless steel, aluminium

A = LEFT/RIGHT  
B = UP/DOWN  
C = IN/OUT

Contact us at  
[info@kleindiek.com](mailto:info@kleindiek.com)  
 or find your local agent at  
[www.kleindiek.com](http://www.kleindiek.com)